

CRYO PRO™ & CRYO PRO™ XL

COLD COMPRESSION
THERAPY



COLD COMPRESSION THERAPY



Cold compression therapy works by constricting the blood vessels, which in turn reduces swelling and inflammation. The Cryo Pro™ and Cryo Pro™ XL control unit features automated compression and easy set-up for patient comfort.

The Cryo Pro™ and Cryo Pro™ XL Unit features nine different modes for varying pressure and time sessions. You can choose from 3 different modes of pulse compression – low, medium and high – to accommodate your cold compression needs. With its dual hoses it allows you to use two wraps simultaneously.

Cold compression therapy can also calm muscle spasms, minimize bruising, and soothe nerve endings, helping encourage quick and effective healing for a variety of muscle, tendon, and ligament injuries.

In addition, cold compression can help reduce a patient's dependence on painkillers following injury.

- Provides pain relief
- Reduces swelling and inflammation
- Prevents or limits bruising
- Slows or stops bleeding

'encourage quick and effective healing for a variety of injuries'

TECHNICAL INFORMATION

- Intermittent cold water flow recovery device
- 6L Reservoir to be filled with water and ice
- Water temperature adjustable between 33.8°F and 50°F (based on the amount of ice used)
- Gradual compression cycle combined with cold diffusion
- Powerful pump: 700L/h and 25W, allowing simultaneous use of two wraps
- Easy to use with only 4 buttons: power on/off, start/stop, and program selection (+/-)
- 9 Programs, ranging from low to high pressure levels, with cycle durations of 20 to 60 minutes and pressure/release times between 8 and 20 seconds
- Quiet operation and LCD screen
- Weight: 8.8 lbs
- 2-Year Warranty

CRYO PRO™ XL UNIT, WRAPS & ACCESSORIES

PRODUCT CODE	DESCRIPTION
DA0302X101	Dr. Aktive Cryo Pro XL Control Unit
DA0302X201	Dr. Aktive Cryo Pro XL 360 Knee Wrap
DA0302X202	Dr. Aktive Cryo Pro XL Shoulder Wrap
DA0302X203	Dr. Aktive Cryo Pro XL Ankle Wrap
DA0302X204	Dr. Aktive Cryo Pro XL Universal Wrap
DA0302X205	Dr. Aktive Cryo Pro XL Full Leg Wrap
DA0302X206	Dr. Aktive Cryo Pro XL Full Long Leg Wrap
DA0302X207	Dr. Aktive Cryo Pro XL Back and Hip Wrap
DA0302X208	Dr. Aktive Cryo Pro XL Extension Strap
DA0302X301	Dr. Aktive Cryo Pro XL Portable Bag
DA0302X303	Dr. Aktive Cryo Pro XL Split Hose
DA03020031	Dr. Aktive Cryo Pro XL Power Bank 2.0

CRYO PRO™ UNIT, WRAPS & ACCESSORIES

PRODUCT CODE	DESCRIPTION
DA03020001	Dr. Aktive Cryo Pro Control Unit
DA03020002	Dr. Aktive Cryo Pro Shoulder Wrap
DA03020003	Dr. Aktive Cryo Pro 360 Knee Wrap
DA03020004	Dr. Aktive Cryo Pro Ankle Wrap
DA03020005	Dr. Aktive Cryo Pro Universal Wrap
DA03020016	Dr. Aktive Cryo Pro Full Leg Wrap
DA03020017	Dr. Aktive Cryo Pro Full Long Leg Wrap
DA03020006	Dr. Aktive Cryo Pro Portable Bag
DA03020007	Dr. Aktive Cryo Pro Power Bank
DA03020008	Dr. Aktive Cryo Pro Battery Case

COLD COMPRESSION THERAPY

CRYO PRO™ RANGE

The Cryo Pro™ and Cryo Pro™ XL are both simple-to-use products, requiring only ice, water, and power (mains or battery). The wraps are simple and easy to apply - one size fits all with the hook and loop fixing allowing the wraps to be adjusted to fit each patient correctly.

The Cryo Pro™ and Cryo Pro™ XL unit have a straightforward digital display with 9 settings to control session time and compression level. The Cryo Pro™ XL Unit has a higher pressure capacity and larger reservoir than the Cryo Pro™ which allows it to support dual hoses and allows you to use two wraps simultaneously. They are both simple to empty and clean, removing the water from the main unit is easy due to the removable lid. The wraps are easy to drain, using the draining gadget.



WRAPS

The Cryo Pro™ range enables you to target specific areas of pain or swelling with our full range of wraps. These include 360 Knee, Shoulder, Ankle, Universal, Full Leg, Full Long Leg and our newest wrap the Back and Hip wrap.

Cryo Pro™ XL Wraps

The Cryo Pro™ XL features enhanced pressure, leading to the development of specially designed wraps that can withstand this increased performance. If you already own wraps from the original Cryo Pro™, you can still use them with the Cryo Pro™ XL. However, for maximum durability and longevity, the tailored Cryo Pro™ XL wraps are the ideal choice.

Full Leg Wrap

Full Leg Wrap Length: 25 Inches (63.5cm)

Full Long Leg Wrap

Full Long Leg Length: 32.3 inches (82cm)

Portable Bag

The Cryo Pro™ Portable Bag comfortably holds the unit and wraps with an adjustable strap, offering extra storage for accessories like the drain gadget, power bank, hose or dual hose, and additional wraps.



360 KNEE




SHOULDER



ANKLE



BACK AND HIP

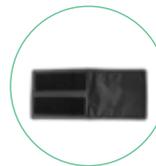

EXTENSION STRAP



FULL LEG



UNIVERSAL

The color of the lining corresponds to the unit its to be used with:

Red lining = Cryo Pro™,

Black lining = Cryo Pro™ XL





Visit us online
www.kinetecusa.com

Inquiries:

Kinetec USA Inc.
N174 W21475 Alcan Drive
Jackson, WI 53037

Tel: (262) 677-1248
Fax: (262) 677-1314
contact@kinetecusa.com
www.kinetecusa.com

dr. aktive[®]
brought to you by
MedFac